


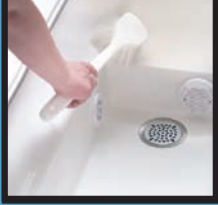




Whirlpool Series (WP6300/WP6000)

ADVANTAGE™ SEATED BATHING SYSTEM








Bath Cleaning Process

Place the chair in the tub, then:

	<p>Turn the black selector knob to "TUB CLEANER."</p> <p>Turn the other black selector knob to the "ON" position.</p>
	<p>Close the drain.</p> <p>Push the button to turn the whirlpool motor on.</p>
	<p>When foamy solution comes out of the jets, wait a moment and turn the black selector knob to "RINSE."</p>
	<p>Use the soapy solution to scrub the tub and chair.</p>
	<p>When the water from the jets runs clear and the tub is completely rinsed, push the button to turn the whirlpool motor off.</p> <p>Turn the black selector knob to the "OFF" position.</p>
	<p>Open the drain.</p> <p>Use the shower wand to rinse the tub and chair.</p>

Bath Disinfecting Process

Place the chair in the tub, then:

	<p>Turn the black selector knob to "DISINFECT."</p> <p>Turn the other black selector knob to the "ON" position.</p>
	<p>Close the drain.</p> <p>Push the button to turn the whirlpool motor on.</p>
	<p>When foamy solution comes out of the jets, allow 2 - 3 inches of liquid to pool in the foot well. Push the button to turn the whirlpool motor off.</p> <p>Turn the black selector knob to the "OFF" position.</p>
	<p>Use the soapy solution to scrub the tub and chair.</p> <p>Leave wet for 10 minutes.</p>
	<p>After 10 minutes, turn the black selector knob to "RINSE."</p> <p>Turn the other black selector knob to the "ON" position.</p> <p>Push the button to turn the whirlpool motor on.</p>
	<p>When the water from the jets runs clear and the tub is completely rinsed, push the button to turn the whirlpool motor off.</p> <p>Turn the black selector knob to the "OFF" position.</p>
	<p>Open the drain.</p> <p>Use the shower wand to rinse the tub and chair.</p>